



# ESTREMO

CITIEFFE NAILING SYSTEM

## Intramedullary nailing system


- ✓ Unique design  
for all lower limb
- ✓ Multiple  
indications
- ✓ Versatile locking  
options
- ✓ Minimal  
instrumentation

**citieffe**<sup>®</sup>  
Essential Moves in Trauma

Surgeons must always decide on the best approach to follow according to their clinical judgement and the patient's needs. Citieffe does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery. The information presented is intended to demonstrate the extent of Citieffe product offerings. Before use, a surgeon must always consult the package insert, product label and/or instructions for use. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Citieffe representative if you have questions about the availability of Citieffe products in your area.

## ESTREMO citieffe nailing system

Intramedullary nailing system with a unique design allowing both an easier tibial suprapatellar approach with semi-extended knee and the application in all lower limb bone segments.

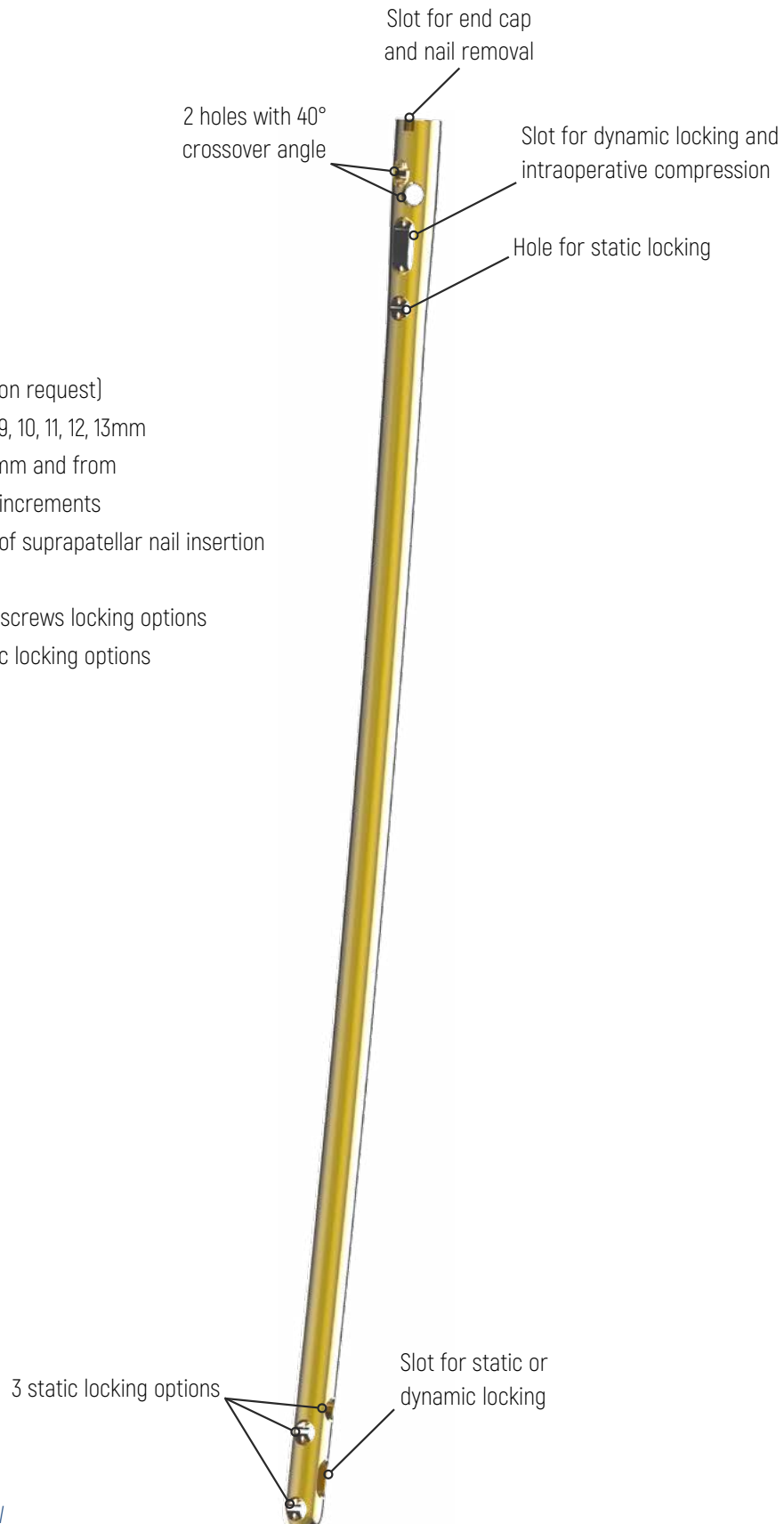

- Titanium alloy ASTM F136
- Solid nail diameter: 8mm (on request)
- Cannulated nail diameter: 9, 10, 11, 12, 13mm
- Available lengths: 150, 200mm and from 260 to 500mm with 15mm increments
- Anatomical bend for ease of suprapatellar nail insertion
- $\varnothing$ 5.2mm screw sizes
- Different multi-directional screws locking options
- Multiple static and dynamic locking options

*Figure shows the cannulated nail*

## Indications

The ESTREMO ctiieffe nailing system is intended for internal fixation of lower limb fractures.

### Femur

The ESTREMO ctiieffe nailing system is indicated for femoral fracture fixation, which may include the following:

- diaphyseal femoral fractures;
- supra-condylar femoral fractures;
- corrective osteotomies;
- pseudoarthrosis;
- non-unions and mal-unions.

### Tibia

The ESTREMO ctiieffe nailing system is indicated for tibial fracture fixation, which may include the following:

- diaphyseal tibial fractures;
- corrective osteotomies;
- pseudoarthrosis;
- non-unions and mal-unions.

### Ankle

The ESTREMO ctiieffe nailing system is indicated for tibiotalocalcaneal arthrodesis.



#### Femur

Antegrade approach  
Retrograde approach



#### Tibia

Suprapatellar approach - Semiextended knee  
Transpatellar approach - Flexed knee



#### Ankle

Tibiotalocalcaneal arthrodesis

---

## ESTREMO instrument set

Only one surgical instrument set with a reduced number of instruments for different implantation techniques.

# ESTREMO

CITIEFFE NAILING SYSTEM

## Intramedullary nailing system



INSPIRED BY PEOPLE  
MOVED BY INNOVATION



Citieffe s.r.l.

Via Armaroli, 21

40012 Calderara di Reno (BO) - Italy

Tel +39 051 721850 - Fax +39 051 721870

info@citieffe.com - www.citieffe.com

YOUR DISTRIBUTOR IS