



**ESTREMO**  
CITIEFFE NAILING SYSTEM

## Intramedullary tibial nail

---

✓ Unique design  
allowing an easier  
suprapatellar approach

---

✓ Multiple options  
to improve  
the stability

---

**citieffe**<sup>®</sup>  
Essential moves in Trauma

Surgeons must always decide on the best approach to follow according to their clinical judgement and the patient's needs. Citieffe does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery. The information presented is intended to demonstrate the extent of Citieffe product offerings. Before use a surgeon must always consult the package insert, product label and/or instructions for use. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Citieffe representative if you have questions about the availability of Citieffe products in your area.

Indicated for the treatment of diaphyseal tibial fractures




- Titanium alloy ASTM F136
- Available lengths: from 260 to 395 mm with 15 mm increments
- Solid nail diameter: 8 mm
- Cannulated nail diameters: 9, 10, 11 mm
- Multiple proximal locking options with the targeting device
- Distal static and dynamic locking options

# Unique design allowing an easier suprapatellar approach

## Semi-extended suprapatellar approach

- — All surgical steps through the protective cannula
- — Simplified nail insertion
- — Guided screw locking

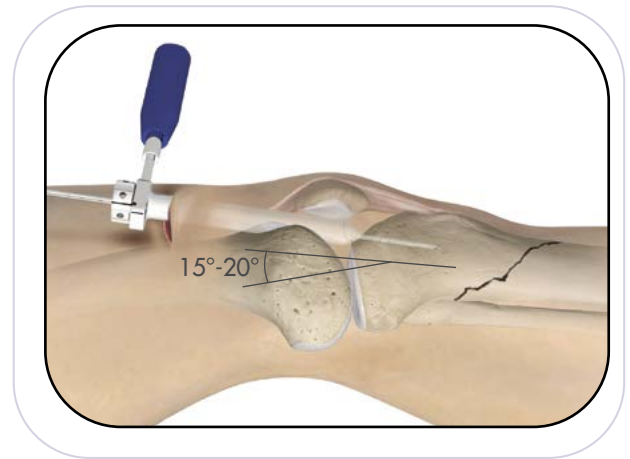

Courtesy of Massimo "Max" Morandi, MD, FACS  
LSU Health Sciences Center - Shreveport, Louisiana USA

# Multiple options to improve the stability

## Multi-directional proximal locking options

- — Two 40° crossover angle screws
- — Two parallel mediolateral screws

## Reduction of bone fragments in tibial plateau

- — Two guided out-nail screws



## Versatile distal locking options

- — Two holes for parallel mediolateral screws
- — Anteroposterior distal hole and slot for static or dynamic locking



Also suitable for transpatellar approach with flexed knee

# ESTREMO

CITIEFFE NAILING SYSTEM

## Intramedullary tibial nail



INTERNAL  
fixation



nails



**citieffe**<sup>®</sup>  
Essential moves in Trauma




Citieffe srl a socio unico

Via Armaroli, 21 - 40012 Calderara di Reno (Bologna - Italy)

Tel +39 051 721850 - Fax +39 051 721870

info@citieffe.com - www.citieffe.com

YOUR DISTRIBUTOR IS



CE  
0476